
Quality standards

Edition: BP 2025 (Ph. Eur. 11.6 update)

Caffeine Hydrate

General Notices

(Caffeine Monohydrate, Ph. Eur. monograph 0268)

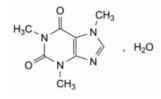

C₈H₁₀N₄O₂,H₂O 212.2 *5743-12-4*

Action and use

Central nervous system stimulant.

Ph Eur

DEFINITION

1,3,7-Trimethyl-3,7-dihydro-1*H*-purine-2,6-dione monohydrate.

Content

98.5 per cent to 101.5 per cent (dried substance).

CHARACTERS

Appearance

White or almost white, crystalline powder or silky crystals.

Solubility

Sparingly soluble in water, freely soluble in boiling water, slightly soluble in ethanol (96 per cent). It dissolves in concentrated solutions of alkali benzoates or salicylates.

It sublimes readily.

IDENTIFICATION

First identification: A, B, E.

https://nhathuocngocanh.com/bp

Second identification: A, C, D, E, F.

- A. Melting point (2.2.14): 234 °C to 239 °C, determined after drying at 100-105 °C.
- B. Infrared absorption spectrophotometry (2.2.24).

Preparation Dry the substance to be examined at 100-105 °C before use.

Comparison caffeine CRS.

- C. To 2 mL of a saturated solution add 0.05 mL of <u>iodinated potassium iodide solution R</u>; the solution remains clear. Add 0.1 mL of <u>dilute hydrochloric acid R</u>; a brown precipitate is formed. Neutralise with <u>dilute sodium hydroxide solution R</u>; the precipitate dissolves.
- D. In a glass-stoppered tube, dissolve about 10 mg in 0.25 mL of a mixture of 0.5 mL of <u>acetylacetone R</u> and 5 mL of <u>dilute sodium hydroxide solution R</u>. Heat in a water-bath at 80 °C for 7 min. Cool and add 0.5 mL of <u>dimethylaminobenzaldehyde solution R2</u>. Heat again in a water-bath at 80 °C for 7 min. Allow to cool and add 10 mL of <u>water R</u>; an intense blue colour develops.
- E. Loss on drying (see Tests).
- F. It gives the reaction of xanthines (2.3.1).

TESTS

Solution S

Dissolve 0.5 g with heating in 30 mL of <u>carbon dioxide-free water R</u> prepared from <u>distilled water R</u>, cool, and dilute to 50 mL with the same solvent.

Appearance of solution

Solution S is clear (2.2.1) and colourless (2.2.2, Method II).

Acidity

To 10 mL of solution S add 0.05 mL of <u>bromothymol blue solution R1</u>; the solution is green or yellow. Not more than 0.2 mL of <u>0.01 M sodium hydroxide</u> is required to change the colour of the indicator to blue.

Related substances

Liquid chromatography (2.2.29).

Test solution Dissolve 0.110 g of the substance to be examined in the mobile phase and dilute to 50.0 mL with the mobile phase. Dilute 1.0 mL of this solution to 10.0 mL with the mobile phase.

Reference solution (a) Dilute 2.0 mL of the test solution to 100.0 mL with the mobile phase. Dilute 1.0 mL of this solution to 10.0 mL with the mobile phase.

Reference solution (b) Dissolve 5 mg of <u>caffeine for system suitability CRS</u> (containing impurities A, C, D and F) in the mobile phase and dilute to 5 mL with the mobile phase. Dilute 2 mL of this solution to 10 mL with the mobile phase.

Column:

- size: I = 0.15 m, $\emptyset = 4.6 \text{ mm}$;
- stationary phase: end-capped amidohexadecylsilyl silica gel for chromatography R (5 μm).

Mobile phase Mix 20 volumes of <u>tetrahydrofuran R</u>, 25 volumes of <u>acetonitrile R</u> and 955 volumes of a solution containing 0.82 g/L of <u>anhydrous sodium acetate R</u> previously adjusted to pH 4.5 with <u>glacial acetic acid R</u>.

Flow rate 1.0 mL/min.

Detection Spectrophotometer at 275 nm.

Injection 10 µL.

Run time 1.5 times the retention time of caffeine.

https://nhathuocngocanh.com/bp

Identification of impurities Use the chromatogram supplied with <u>caffeine for system suitability CRS</u> and the chromatogram obtained with reference solution (b) to identify the peaks due to impurities A, C, D and F.

Relative retention With reference to caffeine (retention time = about 8 min): impurity C = about 0.38; impurity D = about 0.42; impurity F = about 0.6; impurity A = about 0.7.

System suitability Reference solution (b):

— <u>resolution</u>: minimum 2.0 between the peaks due to impurities C and D; minimum 2.5 between the peaks due to impurities F and A.

Limits:

- *unspecified impurities*: for each impurity, not more than 0.5 times the area of the principal peak in the chromatogram obtained with reference solution (a) (0.10 per cent);
- *total*: not more than 0.5 times the area of the principal peak in the chromatogram obtained with reference solution (a) (0.10 per cent);
- *disregard limit*: 0.25 times the area of the principal peak in the chromatogram obtained with reference solution (a) (0.05 per cent).

Sulfates (2.4.13)

Maximum 500 ppm, determined on 15 mL of solution S.

Prepare the standard using a mixture of 7.5 mL of <u>sulfate standard solution (10 ppm SO₄) R</u> and 7.5 mL of <u>distilled</u> water R.

Loss on drying (2.2.32)

5.0 per cent to 9.0 per cent, determined on 1.000 g by drying in an oven at 105 °C for 1 h.

Sulfated ash (2.4.14)

Maximum 0.1 per cent, determined on 1.0 g.

ASSAY

Dissolve 0.170 g, previously dried at 100-105 °C, with heating in 5 mL of <u>anhydrous acetic acid R</u>. Allow to cool, and add 10 mL of <u>acetic anhydride R</u> and 20 mL of <u>toluene R</u>. Titrate with <u>0.1 M perchloric acid</u>, determining the end-point potentiometrically (<u>2.2.20</u>).

1 mL of <u>0.1 M perchloric acid</u> is equivalent to 19.42 mg of C₈H₁₀N₄O₂.

IMPURITIES

Other detectable impurities (the following substances would, if present at a sufficient level, be detected by one or other of the tests in the monograph. They are limited by the general acceptance criterion for other/unspecified impurities and/or by the general monograph <u>Substances for pharmaceutical use (2034)</u>. It is therefore not necessary to identify these impurities for demonstration of compliance. See also <u>5.10</u>. <u>Control of impurities in substances for pharmaceutical use</u>) A, B, C, D, E, F.

A. 1,3-dimethyl-3,7-dihydro-1*H*-purine-2,6-dione (theophylline),

https://nhathuocngocanh.com/bp

B. N-(6-amino-1,3-dimethyl-2,4-dioxo-1,2,3,4-tetrahydropyrimidin-5-yl)formamide,

C. 1,3,9-trimethyl-3,9-dihydro-1*H*-purine-2,6-dione (isocaffeine),

D. 3,7-dimethyl-3,7-dihydro-1*H*-purine-2,6-dione (theobromine),

E. N,1-dimethyl-4-(methylamino)-1H-imidazole-5-carboxamide (caffeidine),

 $F. \quad 1, 7\hbox{-dimethyl-} 3, 7\hbox{-dihydro-} 1 \hbox{\it H-} purine-2, 6\hbox{-dione}.$

Ph Eur